

Millennium™

Lever-Type Amalgam Carrier



Patented Easy-Disassembly Amalgam Carrier!

Improved!

The Millennium lever-type amalgam carrier sets the standard for easy disassembly and fast cleaning of amalgam carriers. No more difficult handling of tiny screws and small screwdrivers to take apart the lever head. With Millennium, you simply push in and hold the Release Button on the side of the handle, gently push down on the curved end until it clears the flat spring, and pull out the end. That's it!

It's an incredibly precise mechanism in an ingeniously clever system, and it makes peaceful coexistence with FDA instrument cleaning recommendations easier than ever before.

Order one today from your favorite full service Dental Dealer, or call Buffalo Dental for more information


Features/Benefits

- **Patented Push-Button Release Heads**
Makes handling, cleaning and complying with FDA recommendations easier than ever before
- **All Stainless Steel Construction**
Helps resist corrosion and increase life
- **Classic Double-Ended Design**
Provides excellent handling and performance

Note: The FDA recommends all surfaces of metal instruments be fully cleaned before autoclaving. This procedure helps ensure the best sterilization result possible. The presence of surface dirt and materials can prevent the steam and elevated temperatures from properly contacting and destroying harmful micro-organisms and bacteria.

Ordering Information

Item #	Description
26810	No. 34 Millennium Break-Apart Amalgam Carrier w/ Nylon Tips
26811	No. 34 Millennium Break-Apart Amalgam Carrier w/ Steel Tips
26815	No. 35 Millennium Break-Apart Amalgam Carrier w/ Nylon Tips
26816	No. 35 Millennium Break-Apart Amalgam Carrier w/ Steel Tips
26820	No. 45 Millennium Break-Apart Amalgam Carrier w/ Nylon Tips
26821	No. 45 Millennium Break-Apart Amalgam Carrier w/ Steel Tips
26811N	No. 34 Regular/Large Millennium Lever-Type Amalgam Carrier w/ TNT coated Tips
26810-3N	"3"- Side Replacement Head Assembly w/ Nylon Tip Complete for Millennium
26810-4N	"4"- Side Replacement Head Assembly w/ Nylon Tip Complete for Millennium
26811-3S	"3"- Side Replacement Head Assembly w/ Steel Tip Complete for Millennium
26811-4S	"4"- Side Replacement Head Assembly w/ Steel Tip Complete for Millennium
88150N	"4"- Side Replacement Head Assembly w/ TNT Tip Complete for Millennium
88155N	"3"- Side Replacement Head Assembly w/ TNT Tip Complete for Millennium

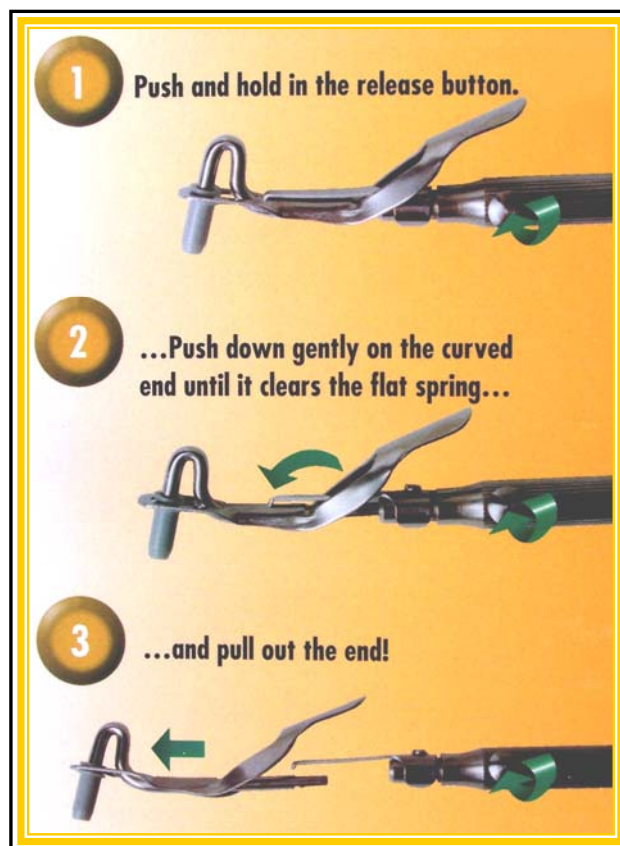

Buffalo Dental Manufacturing Co., Inc.,

159 Lafayette Dr., PO Box 678, Syosset, NY 11791-0678 | Tel: 516-496-7200 or 800-828-0203

Fax: 516-496-7751 | Email: Sales@mailbdm.com | www.buffalodental.com

S/26810 0418